

HOW TO TAKE A SUCCESSFUL INTRAORAL SCAN

Tips from Gnatho

Let Gnathodontics be the digital guide for your office. We believe in creating a cohesive team relationship with your office to create beautiful smiles for your patients.

Full Arch
Scans
Preferred

CLEAR VIEW OF THE MARGIN

The more information we have the better results we can achieve and the less chance of a remake. Make sure the margin is visible with 0.5mm being submarginal with either:

- Double Cord
- Retraction Paste
- Cord Soaked in Hemodent
- Soft Tissue Laser if Needed



MOISTURE CONTROL

- It is important to be proactive with controlling the moisture when taking a digital intraoral scan.
- If space allows and the patient can tolerate it, use cotton rolls and keep the saliva ejector in place as you scan.
- Any liquid that is captured in the scan will either block the margin or reflect the light.

REMOVABLE CASES

Make sure to scan as much information that is possible:

- Full Palate
- All Surfaces of Posterior Teeth
- Capture all Vestibules




IMPLANTS

- Tissue scan is needed with ALL implant cases.
- First, take off the healing cap and scan the tissue immediately.
- Secondly, place the appropriate scan body in place and make sure to scan all sides of the scan body.

SUCCESSFUL CASE

- ✓ Scan while using smooth, gradual movements
- ✓ Retract cheeks, lips and tongue to help capture a successful scan
- ✓ Gnatho will accommodate your office and your preferences. Let us know how to better serve your office and the patient's needs. Preference examples are: occlusion, contacts, cement gap, etc.
- ✓ Double Check Scan
- ✓ Send Scan to cad2cam@gnatho.com

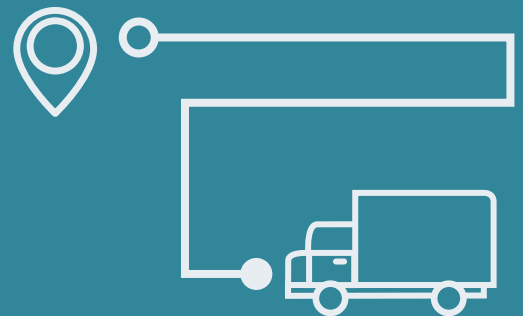

GNATHODONTICS, LTD

— Your Full Service Dental Laboratory —

10488 W. 6th Place | Lakewood, CO 80215

Like • Follow • Share



IOS scanners and how to send to Gnathodontics

Quick Connect: cad2cam@gnatho.com



- Go to Carestream Connect on your scanner
- Search for Gnathodontics or cad2cam@gnatho.com
- Select Gnathodontics when sending cases



- Go to <https://portal.3shapecommunicate.com/login>
 - Search for Gnathodontics or cad2cam@gnatho.com
 - Click Connect and call Gnatho to accept the request
 - Start sending files to Gnatho
- **When sending over a new scan, make a copy file and send



- Call iTero Align Technology at 800-577-8767
- Select option #1
- Use #483 code in order to connect to Gnathodontics
- Select Gnathodontics when sending cases



GNATHODONTICS, LTD

— Your Full Service Dental Laboratory —



- Scan the patient
- Save the files as STL files
- Email the Rx with all of the details needed and the STL files to cad2cam@gnatho.com



- Find a Laboratory on your scanner
- Search for Gnathodontics or cad2cam@gnatho.com
- Click on Request Partnership and call Gnatho to accept the request
- Select Gnathodontics when sending cases

PLANMECA

- Click on My Labs or Find a Lab
- Search for Gnathodontics or cad2cam@gnatho.com
- Make sure the file type is STL
- Use PlanMeca Romexis Cloud to send files to Gnathodontics

Other Intraoral Scanners:

- Scan the patient
- Save the files as STL files
- Email the Rx with all of the details needed and the STL files to cad2cam@gnatho.com

Contact Beth CAD/CAM Manager with any questions: 303.424.9515 (EXT 1750)



GNATHODONTICS, LTD

— Your Full Service Dental Laboratory —

10488 W. 6th Place | Lakewood, CO 80215

Like • Follow • Share

